

MEDICUBE DIGITAL AGE-R



**AGE-R DEVICE
I DEEP SHOT**



**USER
MANUAL**

Quick Start Guide

* When using shot mode and booster mode at the same time, use booster mode after using shot mode.

SHOT MODE



STEP 01

After washing your face, in a dry skin condition rub very lightly along the contour using the head ball.

*Shot mode does not work on wet surface.



STEP 02

Contact the head ball on the skin surface and rub lightly along the contour where you want to improve.

BOOSTER MODE



STEP 01

Apply skincare product on the contour where you want to improve. Turn on the device, and switch the mode to booster mode, then use it.



STEP 02

Gently massage with the head ball on the skin surface, and then apply your daily routine skincare products.

Firm eye contour effect and innovative dual volumizing technology

Use MEDICUBE AGE-R I DEEP SHOT as following.

SHOT MODE



After washing your face, in a dry condition rub very lightly along the contour using the head ball.

BOOSTER MODE



After the SHOT mode, apply skincare product, and massage with the head ball on the skin surface.

Touching Method



SHOT MODE

Lining

Contact the head ball on the skin surface, rub very lightly along the contour where you want to improve.



BOOSTER MODE

Circling

Contact the head ball on the skin surface, massage as if drawing a circle.

Read before use!

It is recommended to start from level 1 when using it for the first time. Please adjust the level appropriately while continuously using it.



Caution

Avoid using around the Adam's apple or on areas with cuts or open wounds of any kind.



Recommended starting level

Start slowly from the lowest setting and find the appropriate level for your skin

SHOT MODE



BOOSTER MODE



I DEEP SHOT Use Cycle

SHOT MODE : 1-2 times a week
BOOSTER MODE : Everyday

Description



Operating Guide

1. ④ Press and hold the power button for 2-3 seconds to turn on the product. Cap *When the power is turned on, the ③ Color LED display lights up with an ON sound.
2. After the product is turned on, it is automatically set to the SHOT mode level 1 with voice guidance, and the mode to be used is set with the ⑤ Mode button. Each time you press the ⑤ Mode button, the mode changes, and each time the mode is changed, the color of the ③ Color LED display changes along with the mode guidance voice.
 - * Shot mode: Red LED | Booster mode: Blue LED
3. After completing the mode setting, use the ④ Level button to adjust the level to be used. ④ Level button, the level goes up one level, and each time the level goes up, the ③ Color LED display flickers along with the level guidance voice.
 - * You can set 1-5 levels, and when using it for the first time, be sure to use it from level 1.
4. After setting the mode and level, touch the ① Head ball to the area of concern.
 - * When using the shot mode, if you press the head ball firmly on the skin, the effect will be reduced, so lightly touch and rub the area of concern without applying force.
 - * Please refer to the detailed usage guide for exact usage.
5. When using the product, the usage time is informed by voice guidance every minute, and it is automatically turned off when the product exceeds 3 minutes. To turn off the device during operation, press the ④ Power button for more than 2 seconds, and the power is turned off with an OFF sound.
 - * When the battery is low, the red light on the ③ Color LED display turns on, and the device is automatically turned off with a low battery sound.
6. For hygienic use, be sure to wipe the ① Head ball with a dry cloth or cotton pad after use, then close the ② Cap and store it.
7. For charging, connect the included C-type charger to the charging port at the bottom of the back of the product.
 - * When charging, the ③ Color LED display turns on, red light is on during charging, and green light is on when fully charged.
 - * The product does not operate while charging.

Contents

Medicube AGE-R I DEEP SHOT,
USB Charging Cable (USB-C),
User Manual



Precautions During Use

1. Dropping the product during use or **lifting** the product by holding only the cap may cause damage and product failure.
2. Make sure no liquid or cosmetic product enters the device.
3. Contact the manufacturer if the product does not work; do not disassemble it.
4. Make sure that the device and the charging cable are connected properly.
5. Before use, make sure that the cables and any accessories are connected properly.
6. Remove any metal accessories that may come into contact with the user's body.
7. Do not use together with electrical appliances such as electric mats.
8. Do not use next to a fire or in hot places, such as under direct sunlight. Keep body and cable away from heat sources.
9. Do not hold the device in one area for long periods of time. Continuously massage in motions.
10. Beware of contact with other products as the device may malfunction.
11. To prevent overheating of internal components, turn off the device for 10 minutes when using the product more than 3 times in a row. Continuous use may shorten life of the device and may cause malfunctions.
12. Do not use it for any purpose other than its intended purpose. It may be the cause of accidents.
13. Please do not use on persons who are paralyzed or unconscious.
14. Do not tap the roller head on solid surfaces, as this may cause malfunction.
15. Do not hit the device against hard objects or drop the device on the floor.
16. Do not use in high temperatures, high humidity, or in areas with excessive dust, such as bathrooms, kitchens, and cars.
17. When using, stay clear of electrical appliances such as other beauty equipment, TV, radio, etc.
18. Do not use the device around the Adam's apple.
19. Do not use on any voltage level other than the voltage level it is rated for. Using above the recommended level may cause excessive heat, fire, short circuit, or malfunction.
20. Do not handle the device and charging cable with wet hands (this may cause electric shock).
21. Do not repair, disassemble or modify the device on your own (this may cause product damage).
22. The device does not operate when charging.

Precautions for User

Please do not use the device if you are in these condition ; particularly bad health condition, drunken condition, taking sleeping pills, feeling severe fatigue and being treated by doctor. If you feel pain or health problem during use, please stop using. If you use the device on same area repeatedly or use in wrong method, it may damage the skin. If you feel something wrong with your skin, keep the affected area cold and consult with specialist. * If your skin is sensitive, the skin can get slightly red even in general use. Please well noted with the manual and precaution and go through appropriate after care.

Do not use if you are in condition down below

1. Anyone whose skin is red and tingling due to excessive tanning.
2. Anyone whose skin is specially sensitive due to allergic constitution.
3. Anyone who is using an electronic device directly worn in body such as cardiac pacemaker.
4. Anyone who is using an electronic medical device for life retaining such as artificial heart-lung machine.
5. Anyone who is using a mounted electronic medical device such as ECG device.
6. Anyone who is using a steroid type hormone drug for long time.
7. Anyone who has telangiectasia due to liver function disorder.
8. Anyone who is outpatient receiving treatment.
9. Area where acne festering and inflamed.
10. Dermatitis area caused by cosmetic product, anyone who has atopic dermatitis.
11. Freckle area accompanying itchiness, freckle caused by physical stimulation.
12. A part of the body that contains metals, plastics, silicon, etc.
13. Hemophilia, heart disease, acute disease, hemorrhagic area, malignant tumor, fracture area, infectious disease, eyelid area, skin area, Children, pregnant women, lactating women and elderly.

* Since the device is not a medical device, do not use for therapeutic purposes.

Precaution During Storage

1. Store at place where it is hard for children to reach and where there is no moisture.
2. Do not store at place with magnetic or electromagnetic wave occurrence. It may cause breakdown.
3. After using the head part, wipe it with wet tissue or clean gauze, soft cloth, etc. And store it in dry condition.
4. Store the body part with the power turned off.
5. Do not wipe the product with inflammable material such as thinner and benzene. It can cause discolor and deformation due to a chemical reaction.

Specification

Model Name	ME-IDST	Product Name	MEDICUBE AGE-R I DEEP SHOT
Manufacturer	Easystem Co., Ltd.	Contact	021-312-1300
Address	Prime Castle Knowledge Industry Center 304, 3F, 25, Podosaon-ro, 118beon-gil, Sinsu-gil, Gyeonggi-do, Republic of Korea		
CE Certificate Registration No.	11-R-PP-146-0317	Method of Use and Precautions	Refer to the user manual
Manufacturing No. / Manufactured Date	Refer to the specification label	Other Necessary Matters	Refer to the user manual

Product Specification

Product Dimensions :

67(L)×11(W)×14(D)mm / 170(L)×28(W)×35(D)mm

Input Voltage 5V1A (USB-C)

Battery

450mAh Lithium-polymer

Power Consumption

maximum of 1W

Weight 120g / 0.27lb

Charge Time about 3 hours

MADE IN KOREA

Charger Specification

Quality Indication in accordance with the Electrical Appliances And Consumer Products Safety Control Act Safety Confirmation Report No. KU-100743-100010

Model Name DTP-102005

Product Name Battery

Rated Voltage 3.7V, 450mAh

Manufacturing Company Shenzhen Data


Power Technology Ltd

Manufactured Country China

Manufactured Date : Indicated separately at the bottom of the product

This device observes Article 15 of the Federal Communications Commission (FCC) regulation and fulfills two conditions below.

- 1) Non-emission of harmful substances designated by the relevant laws.
- 2) Absorbing power of exterior impact that may cause malfunction not intended by the user.

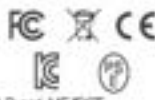

Exchange / Return Policy

1. Warranty

For all Age-R devices, we offer a 1-year limited warranty period from the date of your original purchase. If your device does not work properly due to a manufacturer's defect, we will replace the defective device free of charge within 1 year. This warranty does not cover the devices damaged by: 1) Accident, misuse, or alteration 2) Use with unauthorized accessories or use other than instructed. 3) Connecting to incorrect current and voltage (We only provide C-type USB cable). This warranty only extends if the devices are used in conjunction with authorized device accessories (such as the enclosed charging cables).

2. Returns

For the Age-R devices, all returns must be initiated within 30 days of receiving your order. Please refer to the information as shown below. Simple Remorse: A refund due to a change of mind is only applicable within 30 days of receiving and will only be accepted if the product is returned in original, sealed, and unopened condition. In the case of a return for simple remorse (change of mind), the customer needs to bear the shipping fee. For international orders: You will be responsible for paying for customs fee and returning shipping fee. Additionally, we cannot provide free return shipping labels for international customers.

Damaged or defective product: 1) Please reach out to us via our Customer Service Center (Email) with the following details: First name and last name, email address, phone number, order number. Kindly attach a photo or video describing your issue. 2) We will proceed to issue a confirmation email for you (please note that our Customer Service usually takes 2-3 business days to respond). 3) After receiving a confirmation email, please send the product back to our warehouse via the return address. You will be informed of the return address after the email confirmation. 4) Once the product has arrived at our warehouse, our staff will inspect the condition of the product. 5) In case of product defect, we will proceed to exchange and send a new product to you. *Note: This process may take up to 1 month, including the time taken to ship the product back to our warehouse, and the time taken to inspect the product.

Subject of Recommendation of Lithium Battery Discarding and Charging Method

Warning : Since it may accompany danger such as heating, fire, and explosion, please observe the following matters. 1) Handling Matters Since a swollen battery that can be visually identified may be dangerous, immediately contact the manufacturer or the seller. Use only the designated charger. Keep away from fire (Do not put inside a microwave oven). Do not leave it inside a car in summer. Do not store or use at high temperature and humid place such as sauna. Do not use it for a long time while placing it on blanket, electric pad, and carpet. Do not store at a sealed space while leaving the power on. Be cautious not to contact the electrical terminal with metal products such as needles, coins, key, and watch. 2) Usage The secondary lithium battery warranted by the manufacturer. Do not implement behaviors such as disassembly, compressing, and penetrating. Do not apply abnormal impact such as dropping it from a high place. Do not expose to a temperature above 50 Celsius. Do not contact with humidity. 3) Disposal Guidelines Do not arbitrarily dispose used up battery. When disposing the battery, do not incinerate or put it inside a fire. Deformation and fire may occur due to heat. The method of battery disposal may vary according to countries and regions. Dispose with appropriate method. 3) Recommendations regarding Charging Method Do not charge the battery with the method other than the one shown in the user manual. Fire or explosion may occur. Complete Charge Time : About 3 hours *Charge time may vary depending on the battery remnant amount of portable battery or surrounding temperature or charging environment. It may charge without notice for performance improvement in specifications and exterior appearance.



MEDICUBE DIGITAL AGE-R

So that everyone can enjoy the digital service
which was privilege for few people

